



**THE
BINDING
SITE**

FREELITE™

For Use On
Bayer ADVIA® 1650

Highly Specific Reagents for the Measurement of Immunoglobulin Free Light Chains in Serum

SENSITIVE

Latex enhanced reagents give previously unattainable levels of sensitivity.
<1.0mg/L for serum **free kappa** and **free lambda**¹.

SPECIFIC


Specificity confirmed against whole immunoglobulin preparations. Minimal cross reactivity allows true free light chain levels to be precisely measured in serum despite myelomatous levels of whole immunoglobulin².

SIMPLE

Full application data and detailed instructions on programming are available. All reagents necessary for running the assays, additional to standard system buffers, are supplied in the kits.

RELIABLE

Reference ranges have been established and published in *Clinical Chemistry*³. Many of the world's leading reference laboratories are now routinely offering Freelite assays as part of their Myeloma investigation panel.

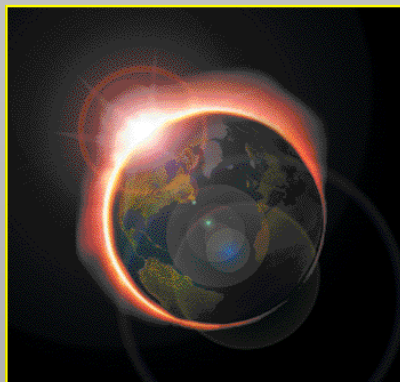

References:

1. Serum Free Light Chain Analysis. AR Bradwell , GP Mead, HD Carr-Smith 2005.
2. Highly Sensitive, Automated Immunoassay for Immunoglobulin Free Light Chains in Serum and Urine. Arthur R. Bradwell, Hugh D. Carr-Smith, Graham P. Mead, Lian X. Tang, Paul J. Showell, Mark T. Drayson and Roger Drew. *Clin. Chem*, 2001;47:4;673-680.
3. Serum Reference Intervals and Diagnostic Ranges for Free κ and Free λ Immunoglobulin Light Chains: Relative Sensitivity for Detection of Monoclonal Light Chains. Jerry A. Katzmann, Raynell J. Clark, Roshini S. Abraham, Sandra Bryant, James F. Lymp, Arthur R. Bradwell and Robert A. Kyle. *Clin. Chem*, 2002;48:9;1437-1444.

Freelite™ is a trademark of The Binding Site Ltd, Birmingham, UK.

Bayer ADVIA® is a registered trademark of Bayer Corporation, New York, USA.

FREELITE™ for use on Bayer ADVIA® 1650

	KAPPA	LAMBDA
NORMAL ADULT SERUM RANGE (95 percentiles)	3.3 - 19.4 mg/L	5.7 - 26.3 mg/L
MEASURING RANGE Standard sample dilution	3.74 - 56.2 mg/L	7.0 - 93.4 mg/L
SENSITIVITY Minimum serum sample dilution (1/1)	0.75 mg/L	0.88 mg/L
STANDARD SERUM SAMPLE DILUTION	1/5	1/8
MINIMUM SERUM SAMPLE DILUTION	1/1	1/1
Sample volume per test (diluted)	5 microlitres	3 microlitres
Minimum sample volume	50 microlitres	50 microlitres
Reagent volume per test	55 microlitres	55 microlitres
Dead volume of reagent per vial	2000 microlitres	2000 microlitres
INTRA ASSAY PRECISION		
Level 1	6.51%	6.28%
Level 2	2.72%	6.42%
Level 3	2.34%	1.81%
INTER ASSAY PRECISION		
Level 1	4.98%	7.07%
Level 2	1.91%	2.90%
Level 3	2.59%	2.96%
OPEN VIAL STABILITY		
Latex reagent - Liquid ready to use	3 months	3 months
Calibrators & Controls - Liquid ready to use	3 months	3 months
KIT CONTENTS		
Latex reagent	2 x 4.0 mL	2 x 4.0 mL
Calibrator	2 x 1.5 mL	2 x 1.5 mL
Low Control	1 x 1.0 mL	1 x 1.0 mL
High Control	1 x 1.0 mL	1 x 1.0 mL
Supplementary Reagent	1 x 25 mL	1 x 25 mL

ANALYSER	DESCRIPTION	PACK	CODE
Bayer ADVIA® 1650	Freelite Kappa Kit	2 x 50 test	LK016.B
	Freelite Lambda Kit	2 x 50 test	LK018.B

All kits are FDA approved for *in vitro* diagnostic use. See product insert for further information. Kits are CE marked for most European countries, please contact us for the latest information.

THE BINDING SITE LIMITED

P.O. Box 11712, Birmingham, B14 4ZB, UK.

Tel: +44 (0)121 436 1000 Fax: +44 (0)121 430 7061 info@bindingsite.co.uk

www.freelite.co.uk www.bindingsite.co.uk

MKG286

PRINTED IN ENGLAND
Nov 05

